

## Product Screener

The EHS Portable Vibratory Screener provides a simple and GMP solution for screening product in a variety of applications. With interchangeable screen sizes and an integrated Ventilation Sleeve to capture dust, this device is a very versatile product for the important task of capturing foreign material.




## Dust Control

The integrated EHS Ventilation Sleeve captures any dust that may be generated by the transfer of dry powders. Additional dust extraction ports can be fitted to capture displaced air in the receiving container. Air flow rates are controlled by manual damper valves.

## Vibratory Screener

Using a simple pneumatic turbine vibrator, the screener allows for quick screening of material during various product transfer processes. Several vibrator

variations and speed settings allow the system to adapt to various applications from V-Shell Blender discharge to Drum to Drum Weighing operations.

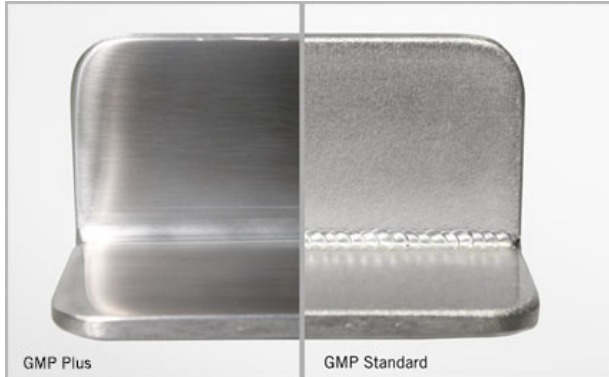


## GMP

This system is fully wash-down capable and can be quickly disassembled for easy cleaning. Complete Stainless Steel construction can be specified in multiple grades from 304 to 316L. Several surface finish options are available as well.

## Easy Screen Change-Out

Multiple screen mesh sizes are available and can easily be switched out. Should the screener not be needed for a particular operation, the screen can be removed and the device can be used as a simple dust containment control.

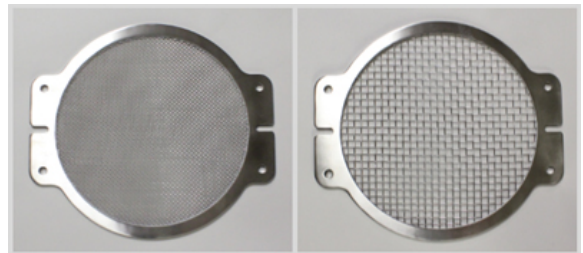

### SHOWS AND EXHIBITS

Interphex 2015  
April 21-23  
Booth #1952

ACHEMA 2015 Germany  
June 15-19

### CONTACT DETAILS

Headquarters:  
1623 W Pioneer Pkwy  
Peoria, IL 61615

Phone: 800-463-7817  
Fax: 309-285-8948

[sales@ehsnow.com](mailto:sales@ehsnow.com)

### LATEST NEWS

EHS now has an application focused sample drawing library for many common industry applications. Visit our Applications page to see some of the applications we can service.

© EHS Solutions LLC

